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"Plasmacluster" and "Device of a cluster of grapes"
are trademark of Sharp Corporation.

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LED irradiation panel
Operation Manual

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Thank you for purchasing this product.
Please read this manual carefully for proper use.
Please be sure to read "SAFETY PRECAUTIONS" before
using the product.
After reading, keep the manual in a convenient location
for future reference.

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Do not use the product for commercial
purposes.

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CONTENTS	Page
FEATURES	34
IMPORTANT SAFETY INSTRUCTIONS	35
PART NAME	38
BASICS OF HOW TO USE	40
· Preparations before use	40
· How to operate	41
MAINTENANCE	44
TROUBLESHOOTING	46
SPECIFICATIONS	47
WARRANTY AND AFTER-SALES SERVICE	48

FEATURES

"Equipped with 4 types of LEDs"

Due to the effects of the lenses that Sharp developed (red LEDs at a wavelength of 630 nm, blue LEDs at 2 wavelengths (415 nm and 455 nm) and near-infrared LEDs at a wavelength of 845 nm), you can focus the light on the front of the panel and use the light efficiently.

What is "Plasmacluster" technology?

Sharp's proprietary technology uses plasma discharge to create and discharge the same positive and negative ions that are in nature. Water molecules that contain Plasmacluster ions adhere to the surface of the skin and form a water molecule coating that retains moisture in the keratin layer and produces the effect of beautiful skin ^{*1}. The ions are the same as those that exist in nature, the mechanism of action has also been clarified ^{*2}, and the safety has been confirmed ^{*3}.

LED x Plasmacluster joint development

"Wake-up mode"

If you set the time you want to wake up, it will wake you up with light and sound. The blue light will gradually brighten from 30 minutes before the set time. The alarm sound can be adjusted to small or large.

Testing institution *1: Sogo Ikagaku Kenkyusho KK
 *2: Tohoku University, Research Institute of Electrical Communication
 *3: LSI Medience Corporation

IMPORTANT SAFETY INSTRUCTIONS

Please comply with this strictly

To prevent injuries or damage to property, it's imperative to observe the following precautions.

■ The consequences of improper use are classified as follows.

 WARNING Indicates the content that may cause death or heavy injury.	 CAUTION Indicates the content that may cause light injury or loss of property.
--	---

■ The rules that should be complied with are indicated with the following "signs" (icons).

 Indicates the prohibited item.	 Indicates to-do.
--	--

WARNING To prevent an accident and an injury by an electric shock and a short circuit, the ignition.

- **Will not be modified. Only qualified personnel should perform disassembly and repair.**
(Doing so could lead to fire, shock, or injury.)

- **Do not dip the product in water or pour water on it.**
- **Do not use this appliance near bathtubs, showers, basins or other vessels containing water.**
(Doing so could lead to shock, short, fire, or failure.)

- **Do not allow a child to use appliance by themselves, Do not store or use within reach of a children.** (Doing so could lead to failure or injury.)
- **Persons who should not use product it should not be used for.**
 - **Persons that are pregnant, Persons that are menstruating**
(may cause skin troubles as hormone balance is unstable)
 - **Persons with head injury, rash, excessive burns, or other skin abnormality and the persons extremely exhausted** • **Persons who are unwell** • **Persons who are drunk** • **Persons extremely exhausted**
 - **Persons who have facial melanosis** • **Persons who have photosensitivity** • **Persons who have epileptic seizures that are triggered by flashing lights** • **Persons who have heart disease**
 - **Persons who have facial neuralgia** • **Persons who have eye problems** • **Persons who have telangiectasis due to steroid hormones or impaired liver function**
- **Persons who should not use areas it should not be used for.**
 - **Areas other than the specified areas** • **Areas where you have had surgery**
- **Do not use with voltage that is not AC 100-240 V.**
- **Do not damage the power cable, and do not use a damaged power cable. Do not modify it, bend it excessively, pull it, or twist it, etc. Do not place heavy objects on it or allow it to be pinched.**
(Doing so could lead to fire, shock.)
- **Do not plug or unplug the AC adapter with wet hands.**
(Doing so could lead to electric shock, short circuits, fire, or injuries)
- **Don't plug the device into multi-outlet extension cords or power strips.**
(Doing so could lead to fire, shock.)
- **Do not use if the AC adapter insert area is loose.**
(Doing so could lead to electric shock, short circuits, or fire)
- **Do not put foreign objects such as pins or wires into the air outlet or air inlet.**
(Doing so could lead to electric shock, injuries, or malfunctions)
- **Use out of the reach of children.**
(Doing so could lead to electric shock, or injury.)

IMPORTANT SAFETY INSTRUCTIONS

Please comply with this strictly

To prevent injuries or damage to property, it's imperative to observe the following precautions.

■ The consequences of improper use are classified as follows.

WARNING Indicates the content that may cause death or heavy injury.

CAUTION Indicates the content that may cause light injury or loss of property.

■ The rules that should be complied with are indicated with the following "signs" (icons).

Indicates the prohibited item.

Indicates to-do.

WARNING To prevent an accident and an injury by an electric shock and a short circuit, the ignition.

- **When the LED program is on, put on the exclusive goggles that were provided, and always close your eyes.** (Otherwise it could lead to eye problems)
When using wake-up mode, you don't need the exclusive goggles that were provided.
- **If you are a new user or if you have sensitive skin, you should limit your use to short time while watching how your skin is affected.** (Otherwise it could lead to skin trouble)
- **Remove glasses and contacts, etc. before use.** (Doing so could lead to accidents or problems)
- **Remove metal items such as earrings or necklaces before use.** (Doing so could lead to accidents or problems)
- **Keep batteries out of the reach of children. If a child swallows anything, consult a doctor immediately.** (Doing so could lead to accidents or problems)
- **Batteries may leak due to improper use, which may cause skin and clothing damage. Observe the following as it may cause stains, rupture or ignition.** (Doing so could lead to accidents or problems)
 - **If leaked battery fluid gets into your eyes, rinse with clean water, and get diagnosed by a doctor quickly to avoid causing impairment.** (Doing so could lead to accidents or problems)
 - **If leaked fluid gets on your body or clothes, rinse it immediately and thoroughly with water.** (Doing so could lead to problems)
 - **Do not place in water or fire, do not dismantle, and do not cause the terminal to short circuit.** (Doing so could lead to accidents)
 - **Do not use rechargeable batteries.** (Doing so could lead to accidents)
- **Be sure to use the included AC adapter** (Doing so could lead to fire, electric shock, and breakdown)
- **Insert the AC adapter plug firmly all the way into the power outlet.** (Doing so could lead to electric shock, smoke, or fire)
- **Regularly remove dust, etc. from the AC adapter.** (Doing so could lead to fire)
- **The following people should always consult with a doctor before using the product.**
 - **Persons with asthma**
(To avoid concerns about symptoms such as coughing, due to the slight amount of air that is continuously blowing)
 - **Persons who are taking medications or feel extremely exhausted**
(Doing so could lead to skin problems or physical problems)
- **If you feel pain in your eyes, skin, or body, or if you feel any abnormalities, stop use immediately and consult a medical specialist.**
(To avoid concerns about eye problems, skin problems, or physical problems)
The condition of your skin may also vary depending on the condition of your body. Please use the product while watching how your skin is affected.

Unplug the product

- **Before doing maintenance, stop operation and always unplug the AC adapter.** (Doing so could lead to electric shock, or injury.)
- **If there are any abnormalities (burning odors, etc.), stop operation and unplug the AC adapter.** (Doing so could lead to fire or electric shock)

CAUTION To prevent an accident and an injury by an electric shock the fire.

- **If you use the batteries in an incorrect way, it may cause breaking or fire. Also, if battery fluid leaks, it may cause the device to corrode and get your hands and clothing, etc. dirty. Make sure to do the following.** (Doing so could lead to failure or problems)
 - Insert the batteries with ⊕ and ⊖ in the correct locations.
 - Do not leave unusable batteries inside the device.
- **If people are near the device when the LED is on, be careful when using it.**
- **For the areas where you have removed hair or shaved, wait a while before using.**

- **Do not use it for any purpose other than skin or scalp maintenance or to wake up (optical alarm).** (Otherwise it could lead to skin trouble on the body.)
- **When the LED is on, don't look directly at the LED light therapy panel.** (Otherwise it could lead to eye problems)
- **When the LED is on, don't put your face too close to the LED light therapy panel.** (Otherwise it could lead to eye problems)
- **Don't drop the main unit or AC adapter, don't hit them against anything, and don't apply strong impacts.** (Doing so could lead to failure or injury.)
- **Don't use the exclusive goggles that were provided unless the LED is on.** (Doing so could lead to injury.)
- **Don't rotate the LED light therapy panel beyond its rotation range.** (Doing so could lead to damage or tipping the product over)
- **Don't put things on top of the main unit or the LED light therapy panel.** (Doing so could lead to failure or injury.)
- **Don't use it in an unstable location. Place it in a solid and stable flat location.** (Doing so could lead to injuries or product damage due to the product falling)
- **Don't use rechargeable batteries (nickel cadmium, nickel metal hydride, etc.).** (Doing so could lead to malfunctions)


Unplug the product

- **When unplugging the AC adapter, always hold the main AC adapter unit while unplugging.** (Doing so could lead to electric shock, short circuits, or fire)
- **If you aren't going to use it for a long time, unplug the AC adapter.** (Doing so could lead to electric shock, electric leakage, or fire due to insulation degradation)

PART NAME

FRONT

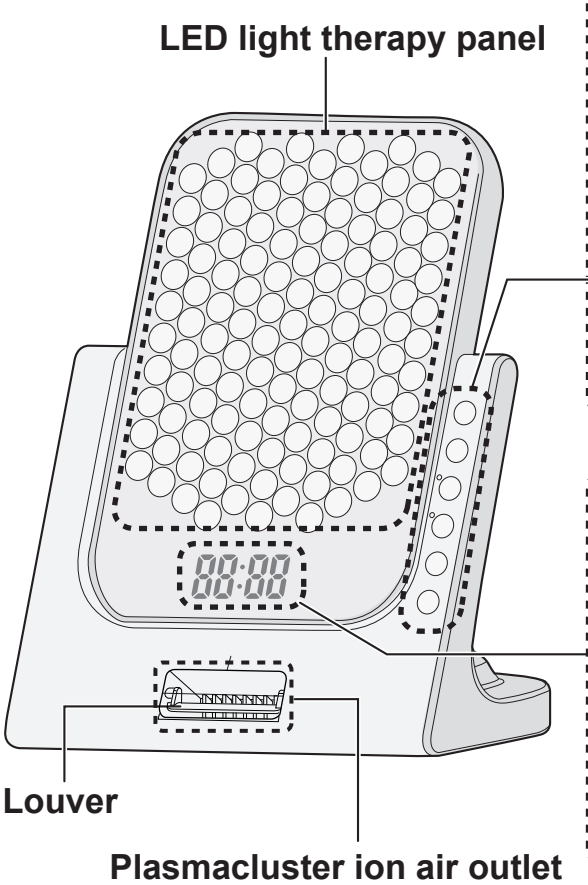LED light therapy panel

Louver

Plasmacluster ion air outlet


Operating area Pages 40 to 43

- LED power button
- Set button
- Clock button [Wake-up light (white)]
- Plasmacluster ion button [Plasmacluster ion light (blue)]
- UP button
- DOWN button

Display area

- ① When setting an LED light therapy program, it displays the program number. Pages 41, 42
- ② When an LED light therapy program is operating, it displays the remaining time that the LED will be on. Page 42
- ③ when setting the time, When setting the time display or when setting wake-up mode, it displays the time. Pages 40, 43

BACK

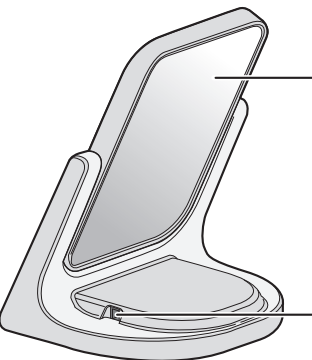
Mirror

If you rotate the LED light therapy panel, you can use it as a tabletop mirror.

Power insertion port


BOTTOM

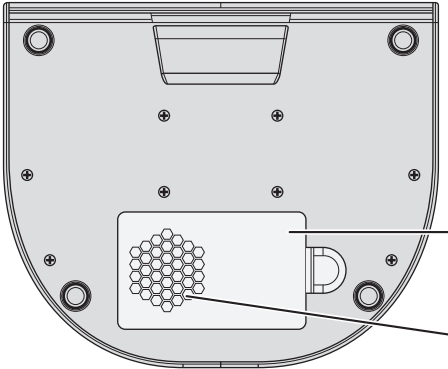Battery cover with filter

Filter part

LED light therapy panel rotation range

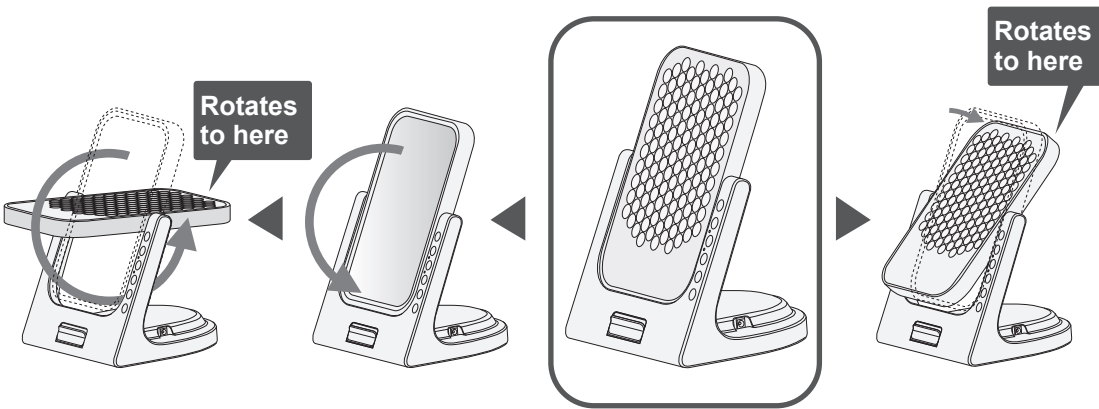Rotates to here


Rotates to here

NOTE

- Do not rotate the LED light therapy panel beyond its rotation range. (Doing so could lead to failure or injury.)
- Do not use force when you rotate the LED light therapy panel. Rotate it slowly. (Doing so could lead to failure or injury.)

ACCESSORIES

- Exclusive goggles
- Unit cleaning brush
- AC adapter
- AAA batteries: 2
- Instruction manual (with warranty)

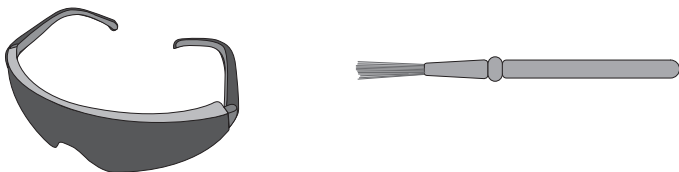
ENGLISH

BASICS OF HOW TO USE

Preparations before use

Insert AAA batteries into the main unit (bottom).

- ① Remove the tab on the battery cover with filter and remove it from the main unit. (Bottom of the main unit)
- ② Insert 2 AAA batteries.
 - Insert the batteries with ⊕ and ⊖ in the correct locations.
 - Insert batteries from the ⊖ side.
- ③ Put the battery cover with filter back on the main unit.

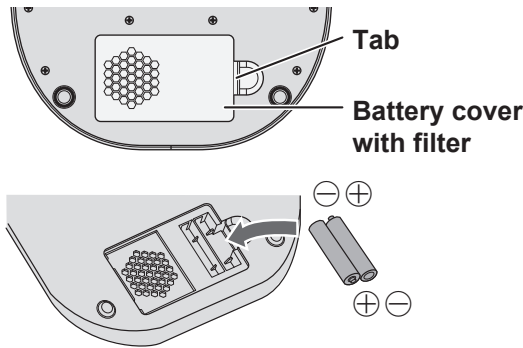
NOTE

- If you don't insert the AAA batteries, the time will not be stored. So if you unplug the AC adapter, the time will be reset.

Connect power to the main unit


- ① Insert the connector into the power insertion port on the main unit.
- ② Insert the AC adapter plug into the power outlet.
 - When pulling out the connector, rotate the connector upward and pull it out.

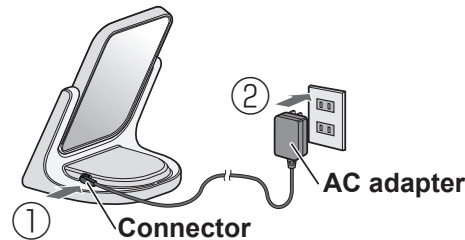Setting the time


- ① If you **hold down** the clock button **for approx. 5 seconds**, the hour display will flash in the display area.

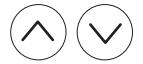- ② Use the up/down buttons to set the hour. If you hold down the up button or down button, it will fast forward through the numbers.

- ③ If you push the clock button, the minute display will flash in the display area.

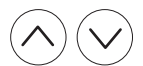
- ④ Use the up/down buttons to set the minute. If you hold down the up button or down button, it will fast forward through the numbers.

- ⑤ If you **hold down** the clock button **for approx. 5 seconds**, the time setting will be complete.

- If the clock cannot be displayed

Please check the TROUBLESHOOTING section. [Page 46](#)

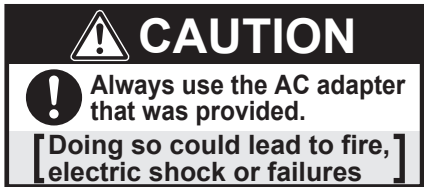How to operate (1)

LED light therapy programs

Program name	LEDs that are lit	Time (minutes)
P-01	Red LED + blue LED + near-infrared LED	10
P-02	Red LED + near-infrared LED	10
P-03	Blue LED + near-infrared LED	10
P-04	Red LED + blue LED	10
P-05	Red LED	10
P-06	Blue LED	10

NOTE

- When the LED is on, put on the exclusive goggles that were provided, and always close your eyes.
- Place it in a solid and stable flat location. (Doing so could lead to injuries or product damage due to the product falling.)


For face

For scalp

Setting silent mode

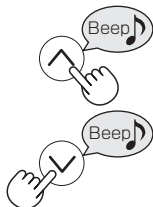You can mute the operating sound by **holding down** the up and down buttons at the same time **for approx. 2 seconds**.

In wake-up mode, the alarm sound will still ring even if silent mode has been set.

- If you **hold down** the up and down buttons at the same time **for approx. 2 seconds** again, the operating sound will turn back on.

BASICS OF HOW TO USE

How to operate (2)

When the LED is on

- ① When you push the LED power button, the program name will be displayed in the display area. 
- ② Push the up/down buttons to select the program.  Page 41
- ③ Adjust the LED light therapy panel angle and the distance from the panel.
Distance from the panel : 10 to 30 cm (For face)
10 to 20 cm (For scalp)
- ④ If you push the set button and count for 5 seconds, the LED will turn on. 
 - While the LED is on, the remaining time will be displayed in the display area. 

<Example>
If the remaining time the light will be on is 7 minutes and 35 seconds.
 - If you push the set button while the LED is on, the LED will temporarily turn off/turn back on. 
 - If you push the LED power button while the LED is on, the LED will turn off. 

• When the LED is on, put on the exclusive goggles that were provided, and always close your eyes.

When using Plasmacluster ions

-  If you push the Plasmacluster ion button, operation will start.
- Each time you push the Plasmacluster ion button, operation will switch in the following order.
- Strong operation**
Plasmacluster ion light turns on

Weak operation
Plasmacluster ion light flashes

OFF
-  Louver
- Please use the louver to adjust the direction the air blows out.

- After you start Plasmacluster ion operation, operation will automatically end after one minute.
- You can use only Plasmacluster ions, or you can use Plasmacluster ions at the same time that the LED is on.
- When the LED is on, put on the exclusive goggles that were provided, and always close your eyes.

NOTE


- If you use cosmetics, etc.* that contain silicone near the main unit, silicone may stick to the inside of the main unit, and it will be more difficult for Plasmacluster ions to be generated. Page 45

* Hair care products (split end coatings, hair mousse, hair treatments, etc.), cosmetics, deodorant, anti-static agents, waterproof sprays, glossing agents, glass cleaners, chemical wipes, waxes, etc.

When using wake-up mode

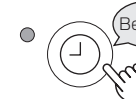[1] Setting the alarm clock

- ① If you push the clock button, the hour display will flash in the display area. While you are setting the alarm clock, the wake-up light will flash.  
 - The default setting is 7:00.
- ② Use the up/down buttons to set the hour. If you hold down the up button or down button, it will fast forward through the numbers. 
- ③ If you push the clock button, the minute display will flash in the display area.  
- ④ Use the up/down buttons to set the minute. If you hold down the up button or down button, it will fast forward through the numbers. 
- ⑤ If you push the clock button, the alarm clock setting will be complete. 
 - If you don't do any operations for approx. 20 seconds, the alarm clock setting will not change, and it will switch to the time display.

[2] Setting wake-up mode

-  If you **hold down** the set button for **approx. 2 seconds**, wake-up mode will be set, and the wake-up light will turn on, and the set time is displayed for approx. 2 seconds. 
- If you **hold down** the set button for **approx. 2 seconds** again, the setting will be cleared.
 - To reduce the volume of the alarm sound
 - If you **hold down** the down button for **approx. 2 seconds**, the reduced volume will be set, and the bottom right of the display area will light up. 

To clear the setting, **hold down** the down button for **approx. 2 seconds**.

• Wake-up mode operation

The LED will turn on 30 minutes before the alarm clock time, and the light will gradually become brighter. When the alarm clock time comes, the alarm sound will ring.

- If 10 minutes pass after the alarm clock time, the LED and the alarm sound will automatically turn off.
- If you push the LED power button or the set button, the LED and alarm sound will turn off.

NOTE

- Please set the time in advance. Page 40
- Do not look directly at the LED when wake-up mode is operating. Also, after waking up with the light and sound, promptly press the LED power button or the set button to turn off the LED.
- If the alarm clock is set to a time that is within 30 minutes of the current time, the LED won't turn on before the alarm clock. The LED will turn on with the alarm clock and the alarm sound will ring.
- Please adjust the LED light therapy panel angle and the distance from the panel. If you place the main unit in a location that is higher than your face, please make sure to place it in a location where it won't fall down onto your face. (Doing so could lead to injuries)

• If you have set wake-up mode but it doesn't operate.

Please check the TROUBLESHOOTING section. Page 46

MAINTENANCE


Routine maintenance

Main unit and accessories

Maintenance timing: When it gets dirty

Wipe the surface with a soft cloth.
If cosmetics, etc. stick to the main unit and you leave them as-is, it may cause discoloration or cracking.

If it is difficult to get these off, wet a cloth with water or a diluted neutral detergent, squeeze it out, and wipe it. After that, wipe it well to make sure that no detergent is left on the unit.

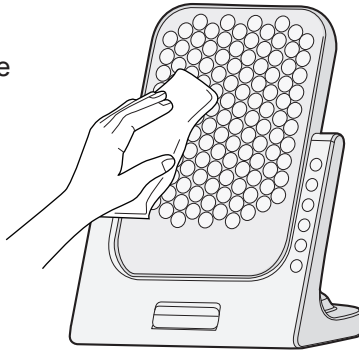Do not use the following products to clean the body and accessories.


(The liquid infiltrating through the discolored area, crack and interior (precision parts) may lead to failure.)

Battery cover with filter


Maintenance timing: About once a week

①Remove the tab on the battery cover with filter and remove it from the main unit.
(Bottom of main unit)

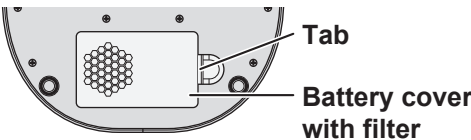②Please use a vacuum cleaner to remove dust from the filter part.

③Put the battery cover with filter back on the main unit.

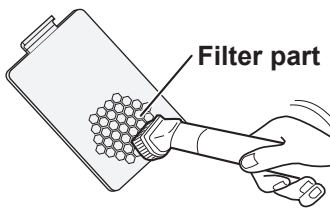NOTE


- Do not push the filter part with strong force.
- The unit cleaning brush that was provided is for cleaning the electrode areas on the unit that generates Plasmacluster ions. Do not use that brush to do filter part maintenance. (Doing so could lead to filter part damage)


Unit that generates Plasmacluster ions

Maintenance timing: About once every 2 months

①Open the louver on the Plasmacluster ion air outlet.

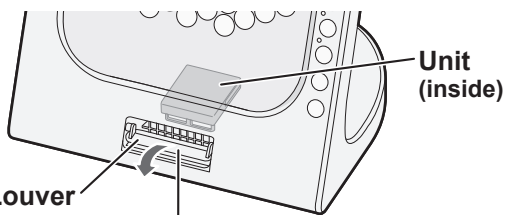②Clean the electrode areas (2 locations) by inserting the unit cleaning brush into the air outlet bar area and moving it forward and backward about 3 times.

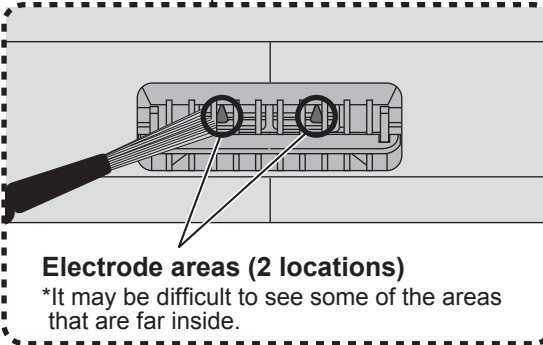
NOTE

- The unit cannot be removed.
- Please do not use anything other than the unit cleaning brush that was provided (do not use cotton swabs, etc.).


Replacing batteries

①Remove the tab on the battery cover with filter and remove it from the main unit.
(Bottom of main unit)

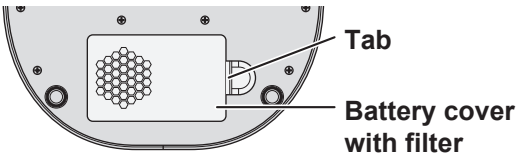②Insert 2 AAA batteries.

- Insert the batteries with ⊕ and ⊖ in the correct locations.
- Insert batteries from the ⊖ side.
- When replacing the batteries, replace both of the old batteries with 2 new AAA batteries of the same type at the same time.

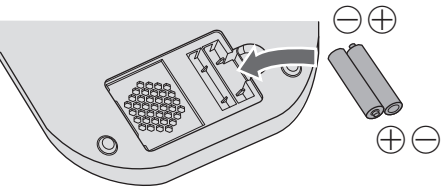③Put the battery cover with filter back on the main unit.

- After you replace the batteries, please set the time again.
- The attached AAA batteries is for checking the operation. The batteries level may be low.
- Don't use rechargeable batteries (nickel cadmium, nickel metal hydride, etc.). (Doing so could lead to operational errors)

TROUBLESHOOTING

Please ask the following questions before commissioning repairs or making consultations.

Problem	Description	Reference
It doesn't operate. I can't get it to work.	<ul style="list-style-type: none">● Is the AC adapter plug connected? ⇒ Please connect it correctly.	40
The time display says 00:00.	<ul style="list-style-type: none">● Are the batteries in? ⇒ Please put the batteries in.● Are your batteries low? ⇒ Please replace them with new batteries, and set the time again.● Did you insert the batteries with the ⊕ and ⊖ in the wrong direction? ⇒ Please check the direction of your batteries.	40, 45
I can't set the time.	<ul style="list-style-type: none">● Isn't LED irradiation in progress? ⇒ Set the time when the LED is not irradiating.● Are your batteries low? ⇒ Please replace them with new batteries, and set the time again.● Did you insert the batteries with the ⊕ and ⊖ in the wrong direction? ⇒ Please check the direction of your batteries.	40, 45
It makes a buzzing sound.	<ul style="list-style-type: none">● That is the sound when Plasmacluster ions are generated, so it's normal.	—
Wake-up mode doesn't work.	<ul style="list-style-type: none">● Has the time been set? ⇒ Please set the time.● Is the alarm clock set to a time that is within 30 minutes of the current time? ⇒ If the alarm clock is set to a time that is within 30 minutes of the current time, the LED won't turn on before the alarm clock. The LED will turn on with the alarm clock and the alarm sound will ring.● Wake-up mode does not operate if a button is pressed or the LED light therapy program is operating at when Wake-up mode is started.	40, 43
There isn't much Plasmacluster ion air coming out.	<ul style="list-style-type: none">● Is the filter part (air inlet) blocked? ⇒ Please remove whatever is blocking it.● Is the air inlet filter part dusty or dirty? ⇒ Please remove the dust or dirt.	44
I am worried about the smell of the Plasmacluster wind coming from the air outlet	<ul style="list-style-type: none">● During the generation of the Plasmacluster ions, ozone are generated, which gives off a kind of smell from the air outlet. The ozone is of the same concentration as natural ozone, so it will not affect health.	—
I can't feel the effect of the Plasmacluster ions.	<ul style="list-style-type: none">● The effects vary for each individual, and the effects may not be immediately evident.● Are the unit's electrode areas dirty? ⇒ Please do maintenance.	45

Displayed error	Description	Reference
E : 01	<ul style="list-style-type: none">● Are you using an AC adapter that is not the one that was provided with this product? ⇒ Please use the AC adapter that was provided.● Please unplug the AC adapter, plug it in again, and try to operate it again. If the error is displayed repeatedly, please stop using it and consult the store.	40
E : 02	<ul style="list-style-type: none">● This error is due to a problem with the Plasmacluster ion blower fan. Please unplug the AC adapter, plug it in again, and try to operate it again. If the error is displayed repeatedly, please stop using it and consult the store.	40

SPECIFICATIONS

Model	AD-BS02, AD-BSC02
Power	100-240V~, 50/60Hz (AC adapter)
Power consumption (When AC adapter is used)	Approx. 25 W
Weight	Approx. 1400 g (excluding the AC adapter, AAA batteries)
Exterior dimensions	(Width) 230 mm (Depth) 194 mm (Height) 303 mm (excluding the AC adapter)
Accessories	AC adapter (length of cable: approx. 1.8 m), exclusive goggles, unit cleaning brush, AAA batteries (2) [AD-BS02 only], instruction manual (with warranty)

WARRANTY AND AFTER-SALES

When you need repairs Take it in for repairs

- 1 Please check the TROUBLESHOOTING section.
Pages 46 to 47
- 2 If you try the suggestions in the TROUBLESHOOTING section and are still having a problem, stop using the product, and always unplug the AC adapter from the power outlet.
- 3 Please ask at the store where you purchased the product.

During the warranty period

- For repairs, please present the warranty. Repairs will be made according to the provisions in the warranty.

If the warranty period has expired

- If you can use the product if it is repaired, if you want it repaired, we will charge you a fee to make the repairs.

Notes for convenience

It's convenient if you write this information here at the time you purchase the product.

Store name

Phone () — (Month) (Day) (Year)

Date of purchase

Structure of repair fees

● Repair fees are composed of technical fees and part costs, etc.	
Technical fees	This fee is for restoring a product that had a failure to normal condition.
Part cost	This fee is the cost of the parts used in the repairs.
Inspection cost	This is a fee for diagnosing a failure. This fee will be charged regardless of whether or not there is a failure.

Have the product you love inspected

If you've been using it for a long time, have the product inspected

Does the product have these symptoms?

- The power cable or AC adapter gets abnormally hot.
- The main unit is deformed, or there is a burning odor.
- There is another problem or failure.

To prevent failures and accidents, stop using the product, unplug the AC adapter from the power outlet, and always ask the store for an inspection. Please consult the store about the costs required for the inspection or repairs.

- Even if the product doesn't have any of the symptoms on the left, if 3-4 years have passed since you made the purchase, for safety reasons, we recommend that you have the product inspected.

SERVICE

Warranty (the warranty is on the back page)

- **The warranty period is one year from the date of purchase (main unit only).**
Fees may be charged even during the warranty period, so please read the warranty thoroughly.
- **Free warranty provision**
 1. If a failure occurs during the warranty period under usage conditions that were in accordance with the precautions in the instruction manual, it will be repaired for free.
 - To request a free repair, please attach a copy of the warranty to the product and ask the store where you made the purchase.
 2. If you move your residence and need information about where to request repairs, please ask the store or agent where you made the purchase.
 3. As a rule, even during the warranty period, fees will be charged for repairs in the following cases.
 - Failures and damage due to mistakes made during use, improper repairs, and modifications.
 - Failures and damage after purchase due to dropping the product, moving the usage location, or transporting the product, etc.
 - Product problems that are due to parts that have a service life (such as the motor), natural wear and tear of consumable parts, wear, or degradation, etc.
 - Failures and damage due to fire, earthquakes, water damage, lightning damage, other natural disasters, pollution, salt damage, gas damage (sulfuric gas, etc.), abnormal voltage, or using voltage (voltage, frequency) that is not the specified voltage, etc.
 - Failures and damage that occurred when the product was left in a vehicle or on a ship, etc.
 - If a copy of the warranty is not attached.
 - If the product serial number, store where the purchase was made, purchase date, and customer's name are not listed on the warranty, or if a receipt that lists this information is not attached.
 - If the text in the warranty has been modified.
 4. The warranty will not be re-issued, so please make sure to store it in a safe place.

* This free warranty provision is a provision that promises free repairs during the period and under the conditions that are expressly stated in the warranty. It does not restrict the customer's legal rights with respect to the party that is issuing this warranty (party that is responsible for the warranty) or other businesses. After the warranty period has passed, if you have questions about repairs, please ask the store where you made the purchase.

Repair notes

ENGLISH

[illegible]

51